


电子显微镜配套设备及耗材    Grids / Support films substrates / Tweezers and small tools / Vacuum coating and evaporation materials / Diamond knives  
 FIB accessories / Filaments and apertures / Consumables kit

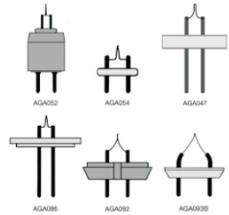

**Micro Hezao**  
 上海禾早电子科技有限公司  
 电话 021-59805057  
 E-mail lfpaull@163.com  
 www.mricohezao.com

Agar Scientific has established an excellent reputation as a leading international supplier of scientific instruments and accessories for over 40 years, specialising in consumables and equipment for all forms of microscopy.

The Agar Scientific servicing team has skilled, certified service engineers who supply, maintain and repair instrumentation and equipment across a range of Electron Microscopy related fields, helping academic institutions and commercial labs across the UK and Europe to keep their equipment operational.




# Agar耗材选型手册

## TEM耗材

<p><b>Standard pattern</b></p> <p>尺寸: 3.05mm</p> <p>材质: 铜、镍、金、铂</p>	<p><b>Tabbed grids</b></p> <p>尺寸: 3.05mm</p> <p>材质: 铜、镍、金</p>
<p>型号: 50mesh, 75mesh, 150mesh, 200mesh, 300mesh, 400mesh</p>	<p>型号: 50mesh, 100mesh, 200mesh, 300mesh, 400mesh, 500mesh, Slot mesh 300*75, Slot2*1mm</p>
<p><b>Hole grids</b></p> <p>尺寸: 3.05mm</p> <p>材质: 铜、镍、金、铂</p>	<p><b>SUPPORT FILMS FORMVAR CARBON</b></p> <p>尺寸: 3.05mm</p> <p>材质: 铜、镍、金</p>
<p>型号: 50µm, 75µm, 100µm, 200µm, 300µm, 400µm, 500µm, 600µm, 800µm, 1000µm, 1500µm, 2000µm</p>	<p>型号: 75mesh, 100mesh, 200mesh, 300mesh, 400mesh, slot2*1mm, F1 finder grid, H7 finder grid</p>
<p><b>Holey carbon films</b></p> <p>孔径尺寸: 约100µm</p> <p>材质: 铜、镍、金</p> <p>型号: 300mesh, 400mesh, H7 finder grid</p>	<p><b>Lacey carbon films</b></p> <p>材质: 铜、镍、金</p> <p>型号: 300mesh, 400mesh, H7 finder grids</p>
<p><b>Square carbon frims</b></p> <p>尺寸: 7*7µm</p> <p>材质: 铜、镍、金</p> <p>型号: 200mesh, 300mesh, 400mesh</p>	<p><b>Circular carbon frims</b></p> <p>尺寸: 1.2µm dia, sparation 1.3µm, 2µm dia, sparation 2µm, 1µm dia, sparation 4µm, 2µm dia, sparation 1wµm</p> <p>材质: 铜、镍、金</p> <p>型号: 200mesh, 300mesh, 400mesh</p>
<p><b>Grid storage boxes</b></p> <p>型号: 50grids, 96grids, 100grids, 50 tabbed grids</p>	<p><b>Cryo grid boxes</b></p> <p>型号: Cryo grid box with lid, round, Cryo grid box base only, round, Cryo grid box with lid, square, Cryo grid box handling rod</p>

## SEM/EDS标样

### SEM多功能校准试样

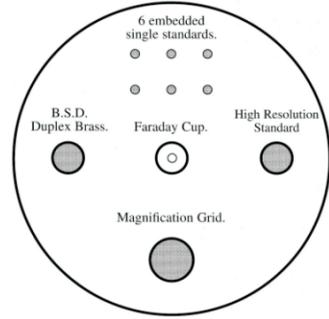

功能:  
6种EDS校准常用标样  
BN, C, Mn, Co, Cu,  
Mn3Al2Si3O12  
B. S. D Duplex brass  
Faraday Cup  
High Resolution Standard  
Magnification Grid




### 高真空仪器标样

仪器类型: SEM, EXCA, Auger  
electron probe


型号: 53种矿物质单质标样

### 金颗粒标样

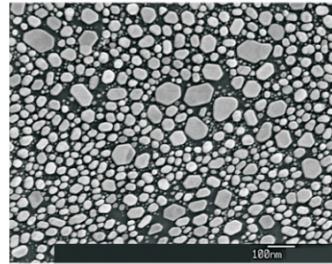

测试范围: 5-150nm

SEM型号:  
JEOL, ZEISS  
TESCAN, FEI  
HITCHE

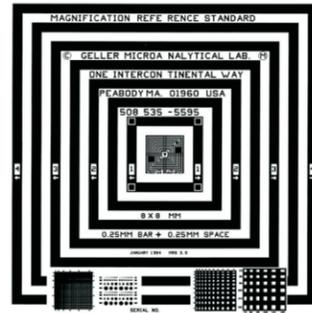

### 标尺校准标样

支持仪器: OM, SEM, TEM  
scanning probe  
LSCM

校准范围: 10X-50000X

## 制样小工具

### 高精度无磁镊子



型号:  
Strong, fine  
flat tips  
very flat tips  
extra fine tips  
Oblique  
straight



### 高精度生物镊子

型号:  
Extra fine, straight tips  
super-fine, very fine  
curved tips

### 高精度自锁镊子



型号:  
fine tips, strong  
flat, rounded tips  
very fine tips,  
curved tips



### 样品台专用镊子

型号: 1/2-3/8  
1/4-3/8  
1/8-3/16  
10mm, 15mm,  
25mm, 32mm

### 真空吸笔

类型: 电池式, 手动

规格: 3mm, 5mm  
7mm, 10mm



### 超薄切片金刚石刀具

型号: 25°、35°、45°  
25nm-200nm  
30nm-1μm  
0.2-5μm

尺寸: 1.5mm, 2mm  
2.5mm, 3mm  
3.5mm, 4mm

## 常用耗材


### 双面碳胶带

型号: 8mm×20m  
12mm×20m  
20mm×20m  
50mm×20m



### 双面碳胶贴

型号: 9mm dia  
12mm dia  
25mm dia

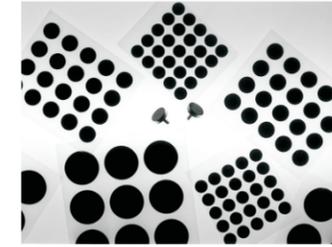

### 导电胶水


规格:  
碳导电胶水  
银导电胶水



### 金属导电胶带

规格: 铝、铜、银

型号:  
6mm×55mm  
12mm×55mm  
25mm×55mm

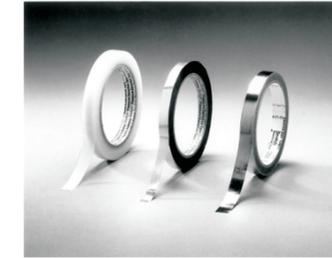

### 进口冷镶嵌树脂

颜色: 无色、黑色  
红色等

特点: 良好收缩比  
无味无毒  
固化时间可控  
4-12H



### 进口冷镶嵌套装

规格: 环氧树脂  
加固剂  
量杯  
模具




### Pon812生物环氧树脂套装

规格: 812树脂  
DDSA硬化剂  
NMA硬化剂  
DMP-30催化剂



### 金刚石砂纸

划痕间隔:  
30μm、15μm  
9μm、6μm  
3μm、1μm  
0.5μm、0.1μm

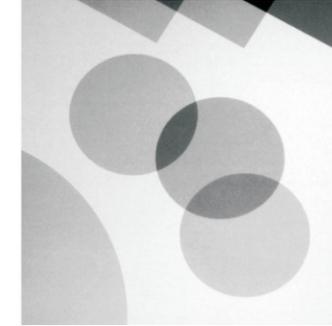

### 抛光砂纸

材质: 碳化硅、氧化铝  
规格: P80、P120  
P240、P320  
P400、P600  
P800、P1200  
P2500、P4000




### 玛瑙研钵

规格: 35×10mm  
50×18mm  
65×18mm



## SEM钨灯丝

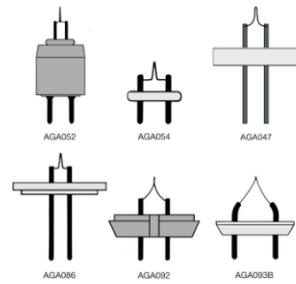

SEM钨灯丝

规格: 预对灯丝  
未对灯丝

型号: TESCAN  
ZEISS  
FEI  
HITCHE  
JEOL



六硼化钨灯丝

规格: 90° cone 15µm microflat  
90° cone 20µm microflat  
60° cone 6µm microflat

型号: TESCAN  
ZEISS  
FEI  
HITCHE  
JEOL

## 离子溅射靶材

碳棒

规格: 6.3mm×150mm  
3.05mm×300mm  
4.6mm×300mm  
6.15mm×300mm



碳绳

规格: 1.55g/m, 2.5mm dia, 3m  
0.7g/m, 1.0mm dia, 1m  
0.4g/m, 1.0mm dia, 1m  
0.27g/m, 0.7mm dia, 5m  
Balzers/Bal-Tec, 3.5m



## SEM/TEM样品台

TEM样品台底座



TEM样品台底座



进口金属靶材

材质: Au, Au/Pa, Pt  
Pt/Pd, Ni, Ag  
Cr等

规格: 57×0.1mm/0.2mm  
54×0.1mm/0.2mm  
58×0.1mm/0.2mm  
40×0.1mm/0.2mm  
115×0.1mm/0.2mm



国产金属靶材

材质: Au, Au/Pa, Pt  
Pt/Pd, Ni, Ag  
Cr等

规格: 57×0.1mm/0.2mm  
54×0.1mm/0.2mm  
58×0.1mm/0.2mm  
40×0.1mm/0.2mm  
115×0.1mm/0.2mm



样品台底座

规格: 可接受定制



钉型样品台

规格: 9.5mm  
45/90° 45°  
70°、45/90°



样品盒

规格: 样品台盒  
高弹膜盒  
钉子样品台盒等



钉型样品台

规格: 12.7mm  
36°、38°  
90°



包埋硅胶板

规格: 10 cavities  
20 cavities  
21 cavities  
24 cavities  
28 cavities



圆弹簧夹片样品台

规格: 平面  
30°、45°  
60°、90°

## FIB配件/耗材

FIB铜网

规格: half-ring  
3 post  
4 post  
5 post



FIB探针

材质: 钨

规格: tip radius < 0.5µm  
10-13° / 6° taper angle



M-Bond 610  
离子减薄胶水套装

规格: resin (14g)  
agent (11g)



三角抛光器

规格: 590 Kid for SEM  
590+ Kit for SEM+TEM  
590S Kit for TEM



抛光膏

规格: 250g



金刚石悬浮液

规格: 0.1µ, 0.25µ  
1.0µ, 3.0µ  
5.0µ, 6.0µ  
9.0µ, 15.0µ  
30.0µ and 45.0µ

# SEM/TEM制样设备

## 电镜相关设备

### 离子溅射仪

Agar auto sputter coater  
Pumping system  
Film thickness monitor  
Chamber size:120\*120mm



### 喷碳仪

Agar auto carbon coater  
Pumping system  
Film thickness monitor  
Chamber size:120\*120mm

### 金刚石慢速切割机

4" (100mm) diameter diamond saw  
6" (150mm) diameter diamond saw



### 亲水化处理仪

Precise and easy vacuum setting  
Short cycle times  
Consistent results  
Auto, programmed and manual mode  
Intuitive touch screen for display and operation  
Supports hydrophilic/hydrophobic and negative/positive modes  
Affordable, compact standalone system  
L-Gasket protects glass jar at base

### 自动抛光仪



### 电解抛光腐蚀仪

电压、电流范围大，可同时满足各种材料的抛光和腐蚀；  
实现恒定电流和恒定电压工作方式；可控制样品的抛光/腐蚀面积，从而保证样品的抛光/腐蚀电流稳定；控制抛光/腐蚀时间；  
搅拌装置保证了抛光/腐蚀介质均匀，样品表面环境一致；  
工作电压、电流可输入计算机，以便于进一步数据分析和研究。

### 冷冻干燥仪

Simple robust construction - easy to maintain.  
Horizontal chamber and large viewing window - excellent visibility of the fluid level and drying process.  
Large robust valves for draining fluids, ingress of CO2 and venting of gas - very durable; the rapid ingress of CO2 helps prevent pre-drying of specimens



### 消磁系统

Magnetic field cancelling systems provide an effective solution for dampening ambient magnetic fields, allowing high resolution performance of TEM, SEM and FEG SEM instruments. The system comprises of a magnetic field control unit, magnetic field sensor and three multicore cables. The cables are installed in the room on the X, Y and Z axes. The sensor measures the strength and phase of the field in real time. The signals are fed to the control unit where three power amplifiers drive anti-phase currents through the multicore cables to effectively nullify the field.

### 精密切割机




### 红外烘烤灯



## EPHEMERON LABS INS 电子束感应电流 ELECTRON BEAM INDUCED CURRENT (EBIC) IMAGING SYSTEMS

We can provide you with a custom solution of hardware and software accessories to make EBIC measurements possible on your SEM. With femtoampere sensitivity, we make it easier to probe current signals that are not accessible with other imaging systems. We do this by providing complete solutions from the most sensitive preamplifier, sample interfacing and shielding, and a purpose designed scan controller with highest resolution possible. All this adds up to measuring very small signals with unprecedented accuracy and very low noise.

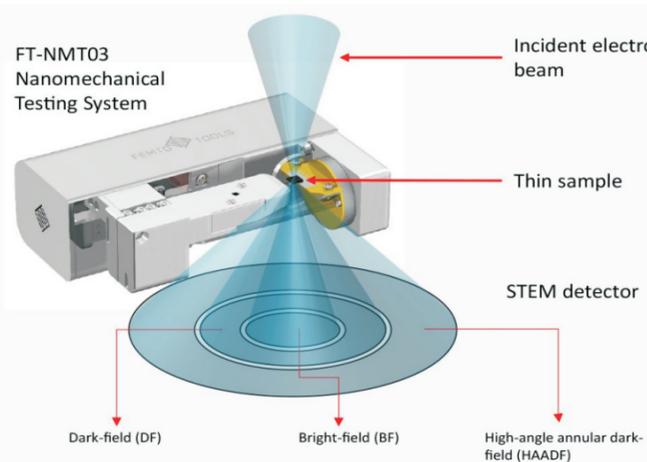

## FEMTO TOOLS

### FT-NMT03 原位纳米力学性能测试系统

能同时实现SEM/FIB高分辨成像和纳米力学性能测试  
力学测量范围0.5nN-200mN(9个数量级)  
位移测量范围0.05nm-21mm(9个数量级)  
能同时实现SEM/FIB高分辨成像和纳米力学性能测试  
能在SEM/FIB最佳工作距离下实现高分辨成像(可达4mm)以及FIB切割和沉积  
五轴(X,Y,Z,旋转,倾斜)位移记录器实现样品台上多样品的自动测试和扫描  
导电的微小传感探针可有效减少荷电效应  
能够通过力和位移两种控制模式实现各种力学测试,例如拉伸、压缩、弯曲、剪切、循环和断裂测试等  
电性能测试模块能够实现力学和电学性能同步测试(样品座配备6个电极)  
实现力学性能测试与其他SEM/FIB原位分析手段联用,如EDX、EBS、离子束沉积和切割  
五轴(X,Y,Z,旋转,倾斜)闭环控制保证样品和微小传感探针的精确对准

## LinDa 离子减薄仪

### SEMPREP2 High-quality site-specific sample preparation in SEM application

focused high-energy ion gun operating from up to 10 keV or optionally ultra high-energy ion gun operating up to 16 keV  
focused low-energy ion gun in the range of 100 eV to 2 keV continuously and independently adjustable milling energy  
beam current density: max. 100mA/cm<sup>2</sup> for focused high-energy ion gun  
max. 150mA/cm<sup>2</sup> for ultra high-energy ion gun  
max. 10mA/cm<sup>2</sup> for focused low-energy ion gun  
sputtering rate: 150 μm/hour on Si at 30° for focused high-energy ion gun  
550 μm/hour on Si at 30° for ultra high-energy ion gun  
28 μm/hour on Si at 30° for focused low-energy ion gun

## ACCURION主动减震台

型号: Sandwich\_1000-4.0  
Sandwich\_1000-6.0  
Sandwich\_1000-8.0

应用领域:  
AFM, STM  
SEM, FE-SEM  
UHV-STM  
纳米力学探针  
共聚焦显微术  
激光扫描显微镜  
光学轮廓仪  
显微硬度检测仪等



WWW.MRICOHEZAO.COM

关于我们  
想知道更多, 请扫描



联系方式: E-mail: lfpaul@163.com TEL:021-59805057